

Genome Assembly Workshop

Deb Triant

University of Virginia

Dept Biochemistry & Molecular Genetics

University of Missouri

Dept Animal Sciences

Programming for Biology

Cold Spring Harbor Labs

17 October 2022



Genome Assembly Workshop

1. Genome assembly with Canu
2. Python problem set



Genome Assembly Workshop

Genome assembly with PacBio data using Canu

Why are we using this program?


- *Relatively* user-friendly and easy to install
- Not computationally intensive with small data sets
- Well maintained and documented
 - version 2.2, released August 2021

Canu assembler

- Canu is a long-read assembler that specializes in PacBio or Oxford Nanopore data
- Derived from Celera assembler
- Corrects reads, trims suspect regions (e.g., adaptors) then assembles the corrected reads
- Can be run with one command to do all steps or each step can be run separately (correcting, trimming and assembling).

A full Canu run includes three stages: correction (green), trimming (red), and assembly (purple).

Generates k-mer histogram and conducts all-vs-all overlaps

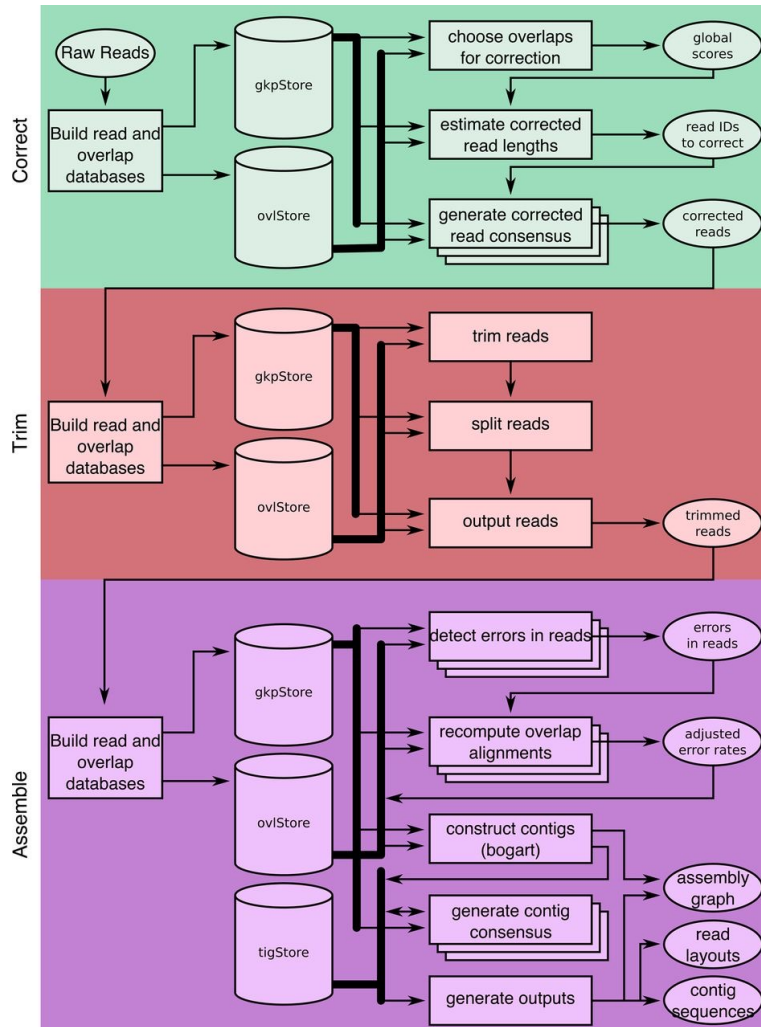

Correct:

Select best overlaps to generate corrected reads. Uses longest 40X reads for correction

Trim:

Identifies unsupported regions and trims reads to longest supported range

Assemble:

Identify sequencing errors, constructs best overlap graph and outputs contigs.

Using Canu

Outputs contigs. No scaffolding.

The unitig construction task finds sets of overlaps that are consistent, and uses those to place reads into a multialignment layout. The layout is then used to generate a consensus sequence for the unitig or a “high-confidence contig”

For eukaryotic genomes, coverage more than 20x sufficient but 30-60x is recommended minimum. Less complex genomes, can work down to 20X coverage.

More coverage, more long reads for Canu to assemble resulting in better assemblies.

Canu Citations

- Koren S, Walenz BP, Berlin K, Miller JR, Phillippy AM. Canu: scalable and accurate long-read assembly via adaptive k-mer weighting and repeat separation Genome Research. (2017).
- Koren S, Rhie A, Walenz BP, Dilthey AT, Bickhart DM, Kingan SB, Hiendleder S, Williams JL, Smith TPL, Phillippy AM. De Novo assembly of haplotype-resolved genomes with trio binning. Nature Biotechnology. (2018).
- Nurk S, Walenz BP, Rhie A, Vollger MR, Logsdon GA, Grothe R, Miga KH, Eichler EE, Phillippy AM, Koren S. HiCanu: accurate assembly of segmental duplications, satellites, and allelic variants from high-fidelity long reads. Genome Research. (2020).

Installing Canu

- We will use the sample Pac Bio *E. coli* data found on the canu homepage: <http://canu.readthedocs.io>
25x *E.coli* fastq files with quality scores.

Found on canu quickstart tutorial:

<https://canu.readthedocs.io/en/latest/quickstart.html#quickstart>

Installing Canu

Getting Help

/canu/documentation/source - help docs

\$ canu command line options

\$ canu -options parameter options

<https://canu.readthedocs.io/en/latest/>

Installing Canu

To install Canu:

Download current release - Canu v2.2

- released August 2021

<https://github.com/marbl/canu/releases>

- click on green “Latest” tab

```
$ wget https://github.com/marbl/canu/releases/download/v2.2/canu-2.2.DArwin-amd64.tar.xz
```

```
#make sure file downloaded 100%
```

```
$ tar -xJf canu-2.2.DArwin-amd64.tar.xz
```

- x extract

- J specific to .xz files

- f file

Installing Canu

Where is it installed?

- Executables found in:
canu-2.2./
**README File
- Add the Canu directory to path:
 - export PATH=/Users/student/canu-2.2/bin:\$PATH
- You can run the assembler with:
 - canu-2.2/bin/canu
- From within bin: ./canu

Canu Sample Data

Download canu test E.coli dataset:

<http://canu.readthedocs.io/en/latest/quick-start.html?highlight=25x>

- 25X subset (223Mb)

Where do you want them to download?

- `wget -O pacbio.fastq`
`https://obj.umiacs.umd.edu/marbl_publications/mhap/raw/ecoli_p6_25x.filtered.fastq`
 - O naming downloaded file `pacbio.fastq`

Canu Sample Data

Check file integrity after downloading

- md5

Confirm the MD5SUM matches:

```
9bb4c10c41c5442d630af8b504042334  pacbio.fastq
```

How many reads within the file?

How many lines?

FASTQ

**@m141013_011508_sherri_c100709962550000001823135904221533_s1_p0/5451
9/26233_32397**

TCGATCGAGTAACTCGCTGCTGTCAGACTGGTTTTTGGTCGATCGACTATTGTTTCAGTCGCAAGAAT
ATTGTGTCCAGTCGATCGACTGAATTCTGCTGTACGGCCACGGCGGATGCACGGTACAGCAGGCTCAG
ACGGATTAAACTGTT

+

5=9=9<=9,-5@<<55>,6+8AC>EE.88AE9CDD7>+7.CC9CD+++5@=-FCCA@EF@+***+*--
55--AA---AA-5A<9C+3+<9)4++=E=+===<D94)00=9))2@624(/(/2/-
(.(6;9((((((.(.'((6-66<6(///

**@m141013_011508_sherri_c100709962550000001823135904221533_s1_p0/5453
2/2817_4395**

GTAAAATTGAGGTAAATTGTGCGGAATTTAGCAATACCGTTTTTTTTTATTATCACCGGATATCTATTC
TGCTGTACGGCCAAGGAGGATGTACGGTACAGCAGGTGCGAACTCACTCCGACGCTCAAGTCAGTGAC
TTAATGATAAGCGTG

+

?????<BBBBBB5<?BFFFFFFECHEFFECCFF?9AAC>7@FHHHHHHHFG?EAFGF@EEDHHDGHHHC
BDFFGDFHF)<CCD@F,+3=CFBDFHBD++??DBDEEEDE:):CBEEEBCE68>?)5?*0?:AE*A
0//:/:*:*:.0)

**@m141013_011508_sherri_c100709962550000001823135904221533_s1_p0/5455
1/25910_41116**

GCTAGTCTTGTGTTTAGTTTTATGTTTTGTCATGTTGTAACGGATTCATAAACATAGGTGTTTGTCTTCT
TTTTATGGTTGTACAATTTGGCCCTAAGGCCCTACACTTACTTGTTTGTCTTTTATGGTACGACAT
TTGAGTGGTGGTTGA

+

Running Canu

- We will run default with correcting, trimming and assembling all in one command:

```
$ canu -p ecoli -d ecoli-pacbio genomeSize=4.8m -pacbio-raw  
pacbio.fastq saveReads=true > CanuRun_20221017.out \  
2>CanuRun_20221017.log **What are we saving?
```

*where are reads....include pathway in your command

-p assembly-prefix

-d output directory

input types:

-pacbio-raw -pacbio-corrected -nanopore-raw -nanopore-corrected fastq or
fasta format

\ = line continuation – don't need if you have room to write it out!!!

****USE THE HELP DOCS!!!!**

Running Canu - Manual assembly

Correct, Trim and Assemble Manually

- You can do the three top-level tasks by hand. This would allow trying multiple construction parameters on the same set of corrected and trimmed reads.

Previous command to correct, trim & assemble:

```
canu -p ecoli -d ecoli-full genomeSize=4.8m -pacbio-raw  
pacbio.fastq saveReads=true 2>full_out_20221017.log **to save  
output
```


Running Canu - Manual assembly

First correct the raw reads:

```
canu -correct -p ecoli -d ecoli-manual genomeSize=4.8m  
-pacbio-raw your_path/canu/pacbio.fastq  
saveReads=true 2>correct_out_20221017.log
```

Then trim the output of the correction:

```
canu -trim -p ecoli -d ecoli-manual genomeSize=4.8m  
-pacbio-corrected ecoli-manual/ecoli.correctedReads.fasta.gz  
saveReads=true 2>trim_out_20221017.log
```

*corrected reads created during first step

Running Canu - Manual assembly

Finally, assemble the output of the trimming twice using two error rates (can use as many as you like):

```
canu -assemble -p ecoli -d ecoli-erate-0.013 genomeSize=4.8m  
correctedErrorRate=0.013 -pacbio-corrected ecoli-  
manual/ecoli.trimmedReads.fasta.gz saveReads=true  
2>assemble-0.013_out_20221017.log
```

```
canu -assemble -p ecoli -d ecoli-erate-0.025 genomeSize=4.8m  
correctedErrorRate=0.025 -pacbio-corrected ecoli-  
manual/ecoli.trimmedReads.fasta.gz saveReads=true  
gnuplotImageformat=svg 2>assemble-0.025_out_20221017.log
```

*trimmed reads created during second step

*The error rate specifies the difference in overlap between two corrected reads which is typically <1% (canu default value 0.045) for PacBio data and <2% (canu default 0.144) for Nanopore data (<1% on newest chemistries). Higher rate, more sensitive.

*Notice in the output there are separate directories for each error rate you specify.

Running Canu

- We will run default with correcting, trimming and assembling all in one command:

```
$ canu -p ecoli -d ecoli-pacbio genomeSize=4.8m -pacbio-raw  
pacbio.fastq saveReads=true > CanuRun_20221017.out \  
2>CanuRun_20221017.log **What are we saving?
```

*where are reads....include pathway in your command

-p assembly-prefix

-d output directory

input types:

-pacbio-raw -pacbio-corrected -nanopore-raw -nanopore-corrected fastq or
fasta format

\ = line continuation – don't need if you have room to write it out!!!

Running Canu

screen - unix command to run programs in multiple windows or “screens”

- If your computer breaks connection, your program will continue run on server

\$ **screen** - to open a new screen

\$ **ctrl-a d** to suspend screen

- Program will continue to run

\$ **screen -r** - to re-open existing screen session

\$ **screen -r #####** - use number to re-open if multiple sessions running

Running Canu - job control

Ctrl-z - stops job running at command line

bg - moves to background

jobs - see what is running

fg - job number - move job to foreground

Running Canu

- We will run default with correcting, trimming and assembling all in one command:

```
$ canu -p ecoli -d ecoli-pacbio genomeSize=4.8m -pacbio-raw  
pacbio.fastq saveReads=true > CanuRun_20221017.out \  
2>CanuRun_20221017.log **What are we saving?
```

*where are reads....include pathway in your command

-p assembly-prefix

-d output directory

input types:

-pacbio-raw -pacbio-corrected -nanopore-raw -nanopore-corrected fastq or
fasta format

\ = line continuation – don't need if you have room to write it out!!!

Canu output files - all found in helpdocs

- **ecoli*.report** - assembly analysis log - histograms of read lengths, kmers and raw/corrected reads. summary of corrected data, overlaps and contig lengths

READS

- **ecoli*.correctedReads.fasta.gz** - The sequences after correction, trimmed and split based on consensus evidence. Typically >99% for PacBio and >98% for Nanopore but it can vary based on your input sequencing quality
- **ecoli*.trimmedReads.fasta.gz** - The sequences after correction and final trimming. The corrected sequences above are overlapped again to identify any missed hairpin adapters or bad sequence that could not be detected in the raw sequences.

SEQUENCE

- **ecoli*.contigs.fasta** - Everything which could be assembled and is part of the primary assembly, including both unique and repetitive elements.
- **ecoli*.unassembled.fasta** - Reads and low-coverage contigs which could not be incorporated into the primary assembly.

How many contigs to we have?

Canu output files - all found in helpdocs

SEQUENCE

- **ecoli.contigs.fasta** - Everything which could be assembled and is part of the primary assembly, including both unique and repetitive elements.
- **ecoli.unassembled.fasta** - Reads and low-coverage contigs which could not be incorporated into the primary assembly.

The header line for each sequence provides some metadata on the sequence:

```
>tig00000001 len=4665109 reads=10220 class=contig suggestRepeat=no  
suggestBubble=no suggestCircular=yes  
trim=20235-4660825
```

len - Length of the sequence **in** bp.

reads - Number of reads used to form the contig.

class - Type of sequence. Unassembled sequences are primarily low-coverage sequences spanned by a single read.

suggestRepeat - If yes, sequence was detected **as** a repeat based on graph topology **or** read overlaps to other sequences.

suggestBubble - If yes, sequence **is** likely a bubble based on potential placement within longer sequences.

suggestCircular - If yes, sequence **is** likely circular.

trim - X-Y to indicate the non-redundant coordinates

Running Canu

Nanopore data

```
canu -p ecoli -d ecoli-oxford genomeSize=4.8m maxInputCoverage=100 -  
nanopore ecolkl2mg1655_R10_3_guppy_345_HAC.fastq
```

PacBio HiFi with HiCanu

```
canu -p asm -d ecoli_hifi genomeSize=4.8m -pacbio-hifi ecoli.fastq
```

Multiple Technologies & Multiple Files

```
canu -p ecoli -d ecoli-mix genomeSize=4.8m -pacbio pacbio.fastq.gz  
-nanopore oxford.fasta.gz
```

Trio Binning Assembly

- Using parental short-read sequencing

```
canu -p asm -d ecoliTrio genomeSize=5m -haplotypeK12 K12.parental.fasta -  
haplotypeO157 O157.parental.fasta -pacbio F1.fasta
```

Genome coverage

- How much coverage is enough?
 - 4.5 Mb - 223Mb fastq reads (25X)
 - 13,124 reads
 - 52,496 lines

Pacific Biosciences released P6-C4 chemistry reads for Escherichia coli K12. You can [download them from their original release](#), but note that you must have the [SMRTpipe software](#) installed to extract the reads as FASTQ. Instead, use a [FASTQ format 25X subset](#) (223MB). Download from the command line with:

Try some assemblies with filtered data sets.

Using canu - split reads

- Count and split files:

```
$ wc -l pacbio.fastq    n=52,496
```

```
$ split -l
```

- I create smaller file with n lines

- 60% 31,498 - lines

```
$ split -l 31498 pacbio.fastq
```

```
$ mv xaa ecoli_filtered_0.60.fastq
```

```
$ rm the rest (xa*)
```

- confirm number of lines and reads

Genome coverage

- How much coverage is enough?
4.5 Mb - 223Mb fastq reads (25X), 13,124 reads
52496 lines

Try some assemblies with filtered data sets:

0.80 - 10,499 reads, 180M, ~22X How many lines?

0.75 - 9,843 reads, 168 Mb, ~18X

0.60 - 7,874 reads, 133Mb, ~15X

0.50 - 6,562 reads, 111Mb, ~12.5X

0.40 - 5250 reads, 89Mb, ~10X

0.25 - 3,281 reads, 56Mb, ~ 6.3X

****Don't forget to split on line counts not reads!**

Running Canu

```
$ canu -p ecoli -d ecoli-pacbio genomeSize=4.8m -  
pacbio-raw pacbio.fastq saveReads=true >  
CanuRun_20221017.out 2>CanuRun_20221017.log
```

****What are we saving?**

- p assembly-prefix
- d output directory

****Make sure you change output file & directory names**

```
canu -p ecoli-60 -d ecoli-pacbio-60 genomeSize=4.8m  
-pacbio-raw ecoli_filtered_0.60.fastq saveReads=true  
> CanuRun-0.6_20221017.out  
2> CanuRun0.6_20221017.log
```

Genome coverage

- How much coverage is enough?

Full data set:

4.8 Mb - 223Mb fastq reads (25X), 13,124 reads

Percent Reads	Coverage	Contigs	Unitigs	Unassembled Sequences	Unassembled Length (Mb)	Genome Size (Mb)
1	25X					
0.8	22X					
0.75	18X					
0.6	15X					
0.5	12.5X					
0.25	6.3X					

Genome coverage

- How much coverage is enough?

Full data set:

4.8 Mb - 223Mb fastq reads (25X), 13,124 reads

Percent Reads	Coverage	Contigs	Unitigs	Unassembled Sequences	Unassembled Length (Mb)	Genome Size (Mb)
1	25X	1	1	840	3.4	4.6
0.8	22X	4	6	639	3.2	4.6
0.75	18X	4	4	566	2.9	4.6
0.6	15X	18	19	641	3.8	4.6
0.5	12.5X					
0.25	6.3X					


Genome coverage

- How much coverage is enough?

Full data set:

4.8 Mb - 223Mb fastq reads (25X), 13,124 reads

Percent Reads	Coverage	Contigs	Unitigs	Unassembled Sequences	Unassembled Length (Mb)	Genome Size (Mb)
1	25X	1	1	840	3.4	4.6
0.8	22X	4	6	639	3.2	4.6
0.75	18X	4	4	566	2.9	4.6
0.6	15X	18	19	641	3.8	4.6
0.5	12.5X	31	35	735	4.6	4.3
0.25	6.3X	67	75	1,188	7.5	1.9

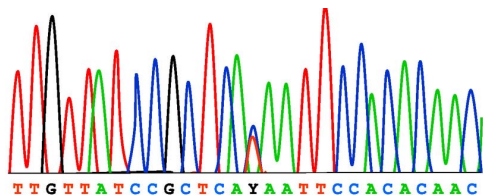

**SANGER
READS**

100%

**HiFi
READS**

Accuracy

80%


**LONG
READS**

0

Read Length (kb)

50

Coverage

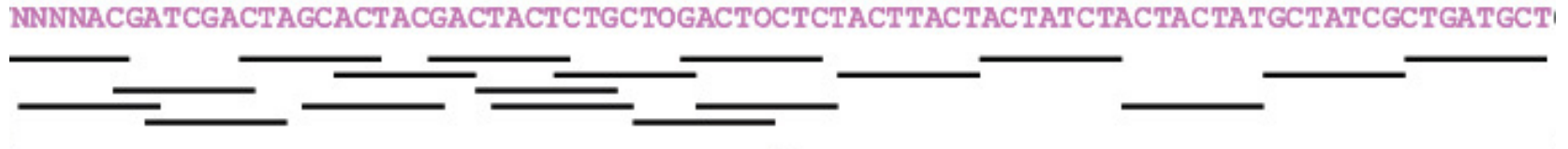

How many times has genome been sequenced?

1. random distribution of reads
2. overlap detection does not vary across genome
3. number of times base sequenced - Poisson distribution
4. can be used to model any discrete occurrence given an average number of occurrences

Poisson Distribution


The probability of a given number of events occurring in a fixed interval of time and/or space if these events occur with a known average rate and independently of the time since the last event.

Reads are not distributed evenly over an entire genome, simply because the reads will sample the genome in a random and independent manner

Lander-Waterman Model

Coverage $c = L * N / G$ where L is the read length, N is the number of reads and G is the genome size

Can be modeled by a Poisson distribution where $\lambda = c$

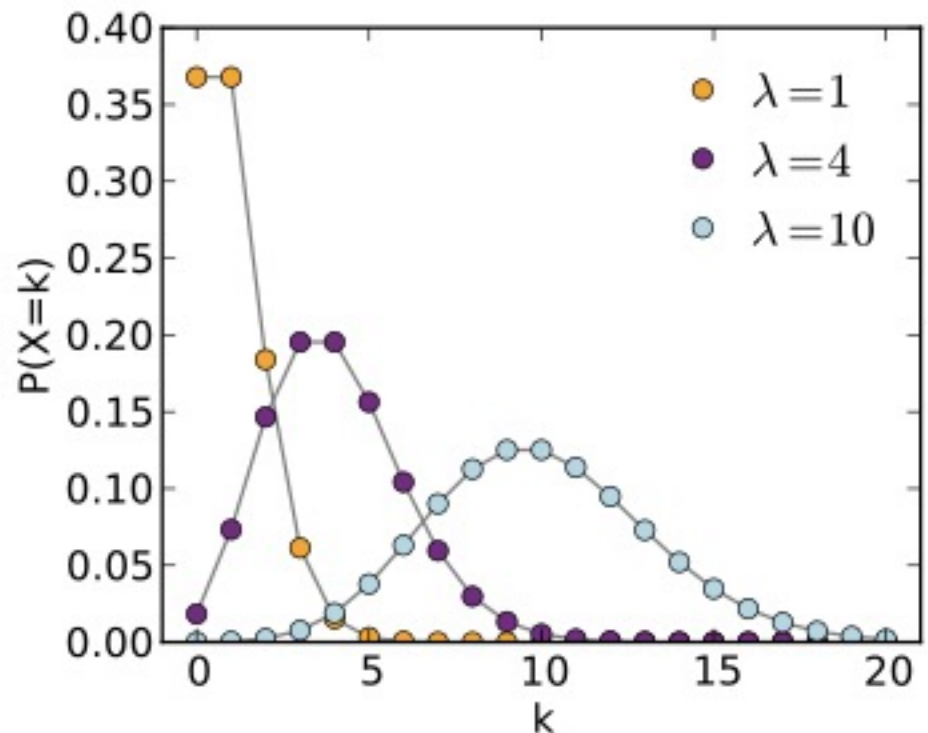

Genome Assembly Problem Set

Using `ecoli-0.25.contigs.fasta`, write a script that reports:

1. The number of contigs in the file
2. The shortest contig
3. The longest contig
4. Total contig length
5. The L50 size
6. The N50 size

The fasta.file for exercises:

[`ecoli_0.25.contigs.fasta`](#).

N50 size

If we place our contigs from largest to smallest on the genome, 50% of the genome in contigs as long as or larger than N50 value

Example: 1 Mbp genome

50%



N50 size = 30 kbp

$(300k + 100k + 45k + 45k + 30k = 520k \geq 500kbp)$

A greater N50 is usually a sign of assembly improvement

- Comparable with genomes of similar size
- Genome composition can bias comparisons
- High L50 vs Low N50

L50 size

Number of contigs that are as long or longer than the N50 value

Example: 1 Mb genome



L50 size = 30 kbp

(300k+100k+45k+45k+30k = 520k \geq 500kbp)

- High L50 vs Low N50
 - longer sequences and fewer of them....in theory
 - lower stringency can inflate N50